



Comfort  
Release®

Sticks when you need it.  
Releases when you don't.



[COMFORTRELEASE.COM](http://COMFORTRELEASE.COM) / 470.636.4410

# ACUTE WOUND CARE/ FIRST AID + ADVANCED WOUND CARE PRODUCTS

DECEMBER 2020 UPDATE

**OUR MISSION** is to improve the clinician and patient experience by providing high quality acute and advanced adhesive wound care products with skin injury prevention solutions.

**OUR VISION** for the lives we touch is to have a compassionate experience with a painless and trauma-free adhesive use and which is easily implemented everywhere care is provided.

Comfort Release® offers a range of **ACUTE WOUND CARE / FIRST AID** and **ADVANCED WOUND CARE** products using the first smart and switchable de-bonding adhesive technology. Our products protect patients from skin pain and skin injuries that commonly occur during adhesive bandage, tape and dressing removal. A simple swipe with rubbing alcohol will release the bond, so you can “turn off” the adhesive bond and remove it without pain, discomfort or injury.

Comfort Release products are **PERFECT FOR CHILDREN, OLDER ADULTS**, anyone taking medications that can affect your skin (such as blood thinners and steroids), or anyone who simply hates the pain of ripping off bandages only to have glue left on the skin.

## ACUTE WOUND CARE - FIRST AID

### COMFORT RELEASE® ADHESIVE BANDAGES & TAPES

*For sensitive skin, easy to apply, water resistant, breathable, long lasting, with painless release. Bandages are made with a nonwoven film and tapes are made with a transparent, bidirectional tear polypropylene film – both with a strong acrylic (latex free) adhesive.*

GB101	1" x 3" Bandage - Single Box w/20 Bandages & 20 Alcohol Prep Pads (20 boxes per case)
GB101-02	1" x 3" Bandage - Single Box w/100 Bandages, NO Alcohol Prep Pads (12 boxes per case)
GB102	2" x 4" Bandage - Single Box w/10 Bandages & 20 Alcohol Prep Pads (20 boxes per case)
GB102-02	2" x 4" Bandage - Single Box w/100 Bandages, NO Alcohol Prep Pads (10 boxes per case)
GB105-02	1" x 10' Tape - Single Blister Pack (50 units per case)
GB105-05	1" x 10' Tape - Box of 12 Rolls (10 boxes per case)
GB106-05	2" x 10' Tape - Box of 6 Rolls (10 boxes per case)
GB107	1" x 3.3' Single Use Tape - Box of 125 Rolls (8 boxes per case)



## ADVANCED WOUND CARE

### COMFORT RELEASE® ADHESIVE TRANSPARENT DRESSINGS

The center transparent windows are made of thin polyurethane material with a low level of adhesion, are waterproof and highly breathable for better wound healing and skin health. The border areas are made of a nonwoven film with a strong acyclic adhesive with our triggered easy release. Made in USA.


GB110	2.375" x 2.75" Transparent Dressing - Box of 10 Single Dressings & 20 Alcohol Prep Pads (20 boxes per case)
GB110-05	2.375" x 2.75" Transparent Dressing - Box of 3 Single Dressings & 6 Alcohol Prep Pads (20 boxes per case)
GB112	4" x 4.75" Transparent Dressing - Box of 10 Single Dressings & 20 Alcohol Prep Pads (20 boxes per case)
GB112-05	4" x 4.75" Transparent Dressing - Box of 3 Single Dressings & 6 Alcohol Prep Pads (20 boxes per case)
GB113	4.5" x 8" Transparent Dressing - Box of 5 Single Dressings NO Alcohol Prep Pads Included (20 boxes / case)



### COMFORT RELEASE® BORDER ADHESIVE FOAM DRESSINGS

Hydrophilic foam pad is made of soft foam, which is comfortable, conformable and highly absorbent. It is breathable and capable of expanding up to 50% when fully saturated. Made in USA.

GB115	4.25" x 4.25" Foam Dressing - Box of 10 Single Dressings & 20 Alcohol Prep Pads (20 boxes per case)
GB115-01	4.25" x 4.25" Foam Dressing - Box of 10 Single Dressings NO Alcohol Prep Pads Included (20 boxes per case)
GB115-05	4.25" x 4.25" Foam Dressing - Box of 3 Single Dressings & 6 Alcohol Prep Pads (20 boxes per case)
GB116	8" x 8" Sacral Foam Dressing - Box of 5 Single Dressings NO Alcohol Prep Pads Included (20 boxes / case)

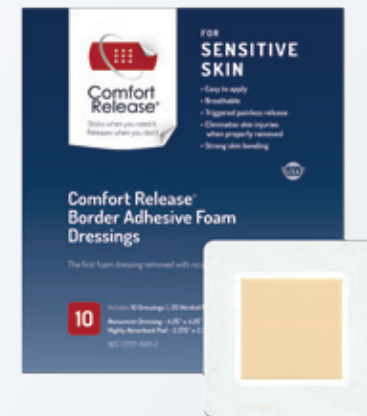

## COMFORT RELEASE® INFUSION KITS

Peripheral and central line infusion kits now with Comfort Release painless, trauma free tapes and transparent dressings are packaged sterile and latex free.

- GB130 IV Kit with 2.375" x 2.75" Comfort Release Transparent Dressing, tourniquet, mini CHG swab stick, Comfort Release Tape roll, 2" x 2" gauze pad, direction card
- 
- GB131 Central Line IV Kit with 4" x 4.75" Comfort Release Transparent Dressing, two 4" x 4" gauze pads, large CHG swab stick, mask, sterile gloves, tissue poly drape, direction card and Guard-IVa CHG dressing pad  
**Expected Release: January 2021**



## WOUND CARE THERAPEUTICS

### COMFORT RELEASE® SILVER ALGINATE DRESSING PADS



The silver alginate dressing is a sterile, single-use dressing pad composed of calcium alginate and silver particles, which absorbs wound exudate and release silver ions in the dressing. As wound exudate is absorbed, the alginate forms a gel, which assists in maintaining a moist environment to help wound healing, and allows intact

removal. Total maximum amount of silver particle is 0.14mg/cm. The silver alginate dressing has broad spectrum antibacterial effects, the main antibacterial effect is silver particles. And the silver protects the dressing from a broad spectrum of bacteria up to seven days.

- GB110-S2 2.375" x 2.75" Transparent Dressing - Box of 10 Single Dressings & 2"x 2" Silver Alginate Dressing Pads - Box of 5 Single Pads & 20 Large Alcohol Prep Pads
- 
- GB112-S4 4" x 4.75" Transparent Dressing - Box of 10 Single Dressings & 4"x 4" Silver Alginate Dressing Pads - Box of 5 Single Pads & 20 Large Alcohol Prep Pads
- 
- GB113-S4 4.5" x 8" Transparent Dressing - Box of 10 Single Dressings & 4"x 4" Silver Alginate Dressing Pads - Box of 5 Single Pads & 20 Large Alcohol Prep Pads

### COMFORT RELEASE® SILVER FOAM DRESSING PADS



The silver foam dressing includes polyurethane foam, polyvinyl alcohol, silver sulfate and silver chloride. It is a highly absorbent foam dressing with 0.35-0.40mg/cm<sup>3</sup> ionic silver as the active component. The dressing provides an optimal moist wound healing environment, combining an effective antibacterial activity in the


dressing with exudate management. The dressing has broad spectrum antibacterial effect for Gram positive and Gram negative bacteria. It has been shown to have 4 log bacterial reduction up to 7 days in vitro against the following bacteria: *Staphylococcus aureus*, *Streptococcus pyogenes*, *Enterococcus faecalis*, *Escherichia coli*, *Klebsiella pneumoniae*, and *Pseudomonas aeruginosa*.

- GB112-F4 4" x 4.75" Transparent Dressing - Box of 10 Single Dressings & 4"x 4" Silver Foam Dressing Pads - Box of 5 Single Pads & 20 Large Alcohol Prep Pads



Use your phone's camera to scan the QR code and watch a video of the removal process.

## MANUKAMED® DRESSING PADS

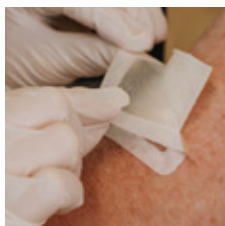

*What is MANUKAMED Honey? Leptospermum honey acts through osmosis and it is thought that its low pH (3.5–4.5) helps modulate the pH of the wound, contributing to an acidic environment conducive to wound healing. ManukaMed honey has methylglyoxal (“MG”) as the leading anti-microbial*

*compound along with a number of phenolics that work in synergy with MG. The activity of ManukaMed honey can pull excess interstitial fluid from the wound and promotes autolytic debridement. ManukaMed honey is available in several forms for your specific wound application.*

- GB102-L2 2” x 4” Bandages - Box of 10 Single Bandages & 2” x 2” ManukaHD Super Lite Dressing Pads - Box of 10 Pads & 20 Large Alcohol Prep Pads

---

- GB110-L2 2.375” x 2.75” Transparent Dressing - Box of 10 Single Dressings & 2” x 2” ManukaHD Super Lite Dressing Pads - Box of 10 Single Pads & 20 Large Alcohol Prep Pads

---

- GB112-L2 4” x 4.75” Transparent Dressing - Box of 10 Single Dressings & 2” x 2” ManukaHD Super Lite Dressing Pads - Box of 10 Single Pads & 20 Large Alcohol Prep Pads

---

- GB112-L4 4” x 4.75” Transparent Dressing - Box of 10 Single Dressings & 4” x 5” ManukaHD Super Lite Dressing Pads - Box of 10 Single Pads & 20 Large Alcohol Prep Pads

---

- GB112-F2 4” x 4.75” Transparent Dressing - Box of 10 Single Dressings & 2” x 2” Manuka Hydrofiber Dressing Pads - Box of 10 Single Pads & 20 Large Alcohol Prep Pads

---

- GB113-L5 4.5” x 8” Transparent Dressing - Box of 10 Single Dressings & 4” x 5” ManukaHD Super Lite Dressing Pads - Box of 10 Single Pads & 20 Large Alcohol Prep Pads

## COMFORT RELEASE® MEDICAL HYDROGEL DRESSING PADS



*Medical Hydrogel Dressing can powerfully absorb water and release moisture. It can effectively absorb wound exudates as well as provide water content to dry wound, maintain moist environment, thus to promote wound healing. The product is semi-transparent, so it is convenient to observe the healing conditions of wound. It will not stick to wound, and no residues left over.*

- GB110-H2 2.375” x 2.75” Transparent Dressing - Box of 10 Single Dressings & 2” x 2” Hydrogel Dressing Pads - Box of 5 Single Pads & 20 Large Alcohol Prep Pads

---

- GB112-H4 4” x 4.75” Transparent Dressing - Box of 10 Single Dressings & 4” x 4” Hydrogel Dressing Pad - Box of 5 Single Pads & 20 Large Alcohol Prep Pads

## BIOPAD™ COLLAGEN DRESSING PADS



*BIOPAD™ is a primary dressing composed of 100% native equine Type I collagen that can accelerate the closure of hard-to-heal wounds. The dressing protects the wound bed from the outer environment, constituting a barrier against exogenous infective agents. It stimulates the formation of new*

*granulation tissue, the proliferation of fibroblasts, and the deposition of new collagen fibers and absorbs wound exudate and can control minor bleeding.*

- GB120-C2 2” x 2” Collagen Dressing Pad - Box of 3 Single Pads

---

- GB121-C4 4” x 4” Collagen Dressing Pad - Box of 1 Single Pad

## ADDITIONAL SUPPLIES

### LARGE ALCOHOL PREP PADS

Saturated with 70% isopropyl alcohol, these latex-free, 2-ply nonwoven pads are ideal for use to trigger our easy release. The pads are also an antiseptic cleanser.

861 Box of 100 Large Alcohol Prep Pads  
(10 boxes per case)



### TECH-MED® ALCOHOL DISPENSING BOTTLES

Manufactured from high density polyethylene. Excellent for dispensing isopropyl alcohol and other cleaners and solvents. Bottles opaque color allows for fluid level visibility. Unique push down dispensing design prevents cross contamination by not allowing dispensed fluid to return to bottle. Locking cap prevents fluid from leaking out.

4024 Single 9oz Alcohol Dispensing Bottle  
(12 bottles per case)



### COTTON BALLS

100% cotton, latex free, non-sterile, soft and uniform in size.

802 Bag of 1000 Cotton Balls  
(2 bags per case)



### STING RELIEF PADS

Saturated with 60% isopropyl alcohol plus 6% soothing benzocaine, these 2-ply nonwoven pads are ideal for use to trigger our easy release. The pads also provide temporary relief of pain and itching associated with minor burns, scrapes, and insect bites.

856 Box of 200 Sting Relief Pads  
(20 boxes per case)



### STICKERS

Roll of stickers to indicate the use of Comfort Release Products. Made in USA. Provided FREE upon request.



## HOW TO APPLY AND REMOVE COMFORT RELEASE® PRODUCTS

We recommend you clean and remove any skin lotion or creams to allow for strong skin bonding of the bandages, tape and dressing.

### COMFORT RELEASE® BANDAGES

To remove a dry bandage, just swipe with rubbing alcohol – we recommend using one or two alcohol prep pads. We recommend that our bandages be replaced at least once every 3 days.

### COMFORT RELEASE® TAPE


Comfort Release® Tape is bidirectional tear, transparent and can be easily sized or cut into strips before using. To remove, use several alcohol prep pads or a cotton ball soaked with rubbing/isopropyl alcohol and swipe back and forth across the outer surface of the tape. Tape will lift from the skin surface in about 8-10 seconds depending on the length of the tape.

### COMFORT RELEASE® DRESSINGS

To apply, peel the release paper then apply the dressing. Gently press down on the dressing for strong securement. The transparent dressing can be left in place for up to 7 days and be used as a secondary dressing cover over a primary therapeutic pad. To *Lift and Look* at a wound or a primary dressing pad, simply swipe one side (or one-half) of the white border with several alcohol prep pads. The Comfort Release® Transparent Dressing may be partially lifted to view a wound and any primary therapeutic dressing pad can be examined and/or replaced.

Once completed, allow the rubbing alcohol to evaporate and replace the Comfort Release® dressing on the skin. It will re-adhere.

To Fully Remove: swipe the entire dry white border area with rubbing alcohol – we recommend using at least two alcohol prep pads. If needed add more alcohol for painless release. The clear window film can be removed without alcohol application and with no discomfort due to its low adhesive skin bond.

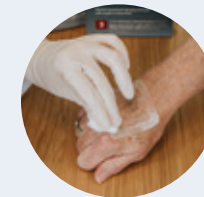

**WARNING:** If irritation or redness develops or persists, discontinue use and consult a physician. People with prior skin reactions to acrylic adhesives, skin burns, or allergic reactions to rubbing alcohol should not use Comfort Release.



Dedicated to Patsy Carter-Rattigan — our #1 patient who graciously volunteered to test Comfort Release samples during our development.

---

Global Biomedical is a fast-growing medical technology company with innovations in chemistry and biomaterials. The company has invested extensively in research and development to find innovative solutions to today's problems and to redefine existing markets.

We developed the Comfort Release® technology to help people of all ages who have experienced the discomfort or injury associated with the removal of bandages, medical tapes or dressings.

Comfort Release is recommended by doctors and nurses, caregivers and other patients - particularly to pediatric, geriatric and any patient with sensitive skin.

To order **ACUTE WOUND CARE / FIRST AID**  
and **ADVANCED WOUND CARE** products  
with Comfort Release technology, call **470.636.4410**.



A Global Biomedical Technology Company

12870 Trade Way Four / #107-314, Bonita Springs, FL, 34135

**COMFORTRELEASE.COM / 470.636.4410**